

SPECIFICATIONS

Lenses: Frames are shipped with clear plastic lenses (lenses are not optical grade). Prescription lenses must be purchased from your eye care professional.

Accessories: Frames ship with one single-lumen and one double-lumen connecting tube as well as twelve nasal prongs (six pairs). With double-lumen tubing user may experience oxygen flow in only one J tube. This is normal and is a function of the oxygen supply or oxygen regulator.

Flow Rate: The maximum flow rate is 5 liters per minute. When using extension tubing - flow rates may vary. Always check with your physician or respiratory therapist to ensure proper oxygenation. The purchase of a finger pulse oximeter should be considered. Check with your physician.

Material: Frames: Metal Tubing and Nasal Prongs: PVC

Weight: Approx. 1.0 oz. (without lenses)


Warranty: One year from date of purchase (frames only).

Colors may vary from those depicted in this booklet. Available styles and colors subject to change.



Nasal prongs and oxygen supply tubing are available in single or double lumen tubing style.

Single lumen used on oxygen systems with one connector. Double lumen used on systems with two connectors.

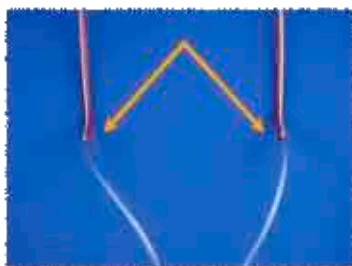

Supply tubing connects to each temple, just aft of the ear pads.



Replaceable nasal prongs (J tubes) connect to the frames near the nose pads.